


Y-Knot® Flex All-Suture Anchor
WITH 2 HI-FI® RIBBON SUTURES

1.3mm Hi-Fi® Ribbon

Double the width and 14% stronger than USP #2 High-Strength Suture¹, CONMED's NEW 1.3mm Hi-Fi® Ribbon combines stronger knot strength with smaller knot stacks. Hi-Fi® Ribbon has a wider format that is designed for broader compression, yet still allows surgeons to tie sliding knots.

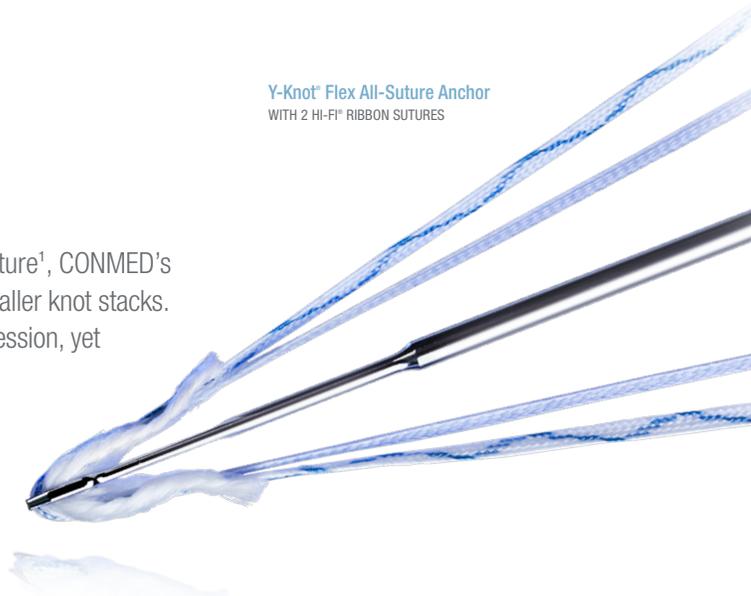

Features and Benefits

Compression

Wide-Format Designed for Broader Compression: With a 1.3mm width, Hi-Fi® Ribbon is twice as wide as #2 Hi-Fi® suture and is designed to provide broader compression for more tendon to bone contact.

Low-Profile

Smaller Knot Stacks:² Designed to help reduce the abrasion and post-operative pain that could be caused by large knot stacks, Hi-Fi® Ribbon provides smaller knot stacks than #2 suture.²

Strength

14% Stronger than #2 High-Strength Suture:¹ In addition to the wider width, Hi-Fi® Ribbon is 14% stronger than #2 High-Strength Suture.¹

Easy-to-Tie

Wide Like Tape, Ties Like #2 Suture: Hi-Fi® Ribbon is designed to provide the width of a tape-style suture while retaining the ability to tie like a #2 suture.

Comfort

Designed to Be Easier on Surgeon's Hands: Compared to the round construction of #2 suture, Hi-Fi® Ribbon is flat and is designed to provide more comfort in surgeons' hands when tying knots.

HI-FI® RIBBON:

- Available pre-loaded on CONMED's Y-Knot® Flex and Y-Knot® RC anchors.



Y-Knot® RC Self-Punching All-Suture Anchor
WITH 3 HI-FI® RIBBON SUTURES

Y-Knot Try It? Schedule a demonstration with your local CONMED Sales Representative or, visit CONMED.com/Ribbon.

¹ Data on File; TR18-00183-1, TR12-441. ² Data on File; TR17-01219.

Ordering Information

	Description	Size	Catalog Number
NEW	Y-Knot® Flex All-Suture Anchor		
	Y-Knot® Flex w/ one 1.3mm Hi-Fi® Ribbon (White/Blue)	1.3mm	Y1301RW
	Y-Knot® Flex w/ one 1.3mm Hi-Fi® Ribbon (Blue)	1.3mm	Y1301RW
	Y-Knot® Flex w/ two 1.3mm Hi-Fi® Ribbons (White/Blue and Blue)	1.8mm	Y1802R
NEW	Y-Knot® RC Self-Punching All-Suture Anchor		
	Y-Knot® RC w/ two 1.3mm Hi-Fi® Ribbon (White/Blue and Blue)	2.8mm	YRC02R
	Y-Knot® RC w/ three 1.3mm Hi-Fi® Ribbon (White/Blue, Blue, and White/Black)	2.8mm	YRC03R