



Hall® MicroFree®

Cordless Small Bone Power System

Unleash Your Precision

The MicroFree® system is designed to make small bone procedures easier. Enhanced ergonomics increase comfort. The freedom of motion optimizes precision. This is small bone power at its finest.




HALL®
POWERED INSTRUMENTS

HALL® TITAN®
PRIMECUT CASSETTE



TRUSHOT® with Y-KNOT®
ALL-IN-ONE SOFT TISSUE
FIXATION SYSTEM

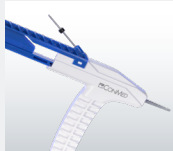

ERGO™
SHOWER HANDPIECE



To learn more about these and other innovative products, call **1-866-4CONMED** or visit **CONMED.com**.

Unleash
Your
Precision.



An Easier Day In The OR.

When we asked surgeons how we could make their jobs easier during small bone cases, the answer was clear – Cordless power handpieces that provide freedom of motion without sacrificing power. These handpieces must be comfortable, ergonomic, well-balanced, and deliver the same power you are accustomed to with corded systems.

Surgeons wanted a new system that delivered improvements without drawbacks and convenience without trade-offs. That's why we developed the MicroFree® Cordless Small Bone Power System.

It's time to cut the cord on your small bone power.

HALL®
MICROFREE®
SMALL BONE
POWER SYSTEM



MicroFree® Pencil Grip



Engineered for
optimal balance

HALL®
MICROFREE®
SMALL BONE
POWER SYSTEM

Cordless Precision.

MicroFree® went through extensive design testing with leading orthopedic surgeons¹ to create a cordless handpiece that makes it easier to perform with precision. Using their feedback, the center of gravity was adjusted to improve balance and control while reducing muscle fatigue.¹ The result was a cordless system that provides freedom of motion in the familiar pencil-grip design.

A Full Cordless Small Bone System.

One Complete System for Your Clinical Needs

The MicroFree® System can be used for most small bone procedures of the foot, ankle, hand, and wrist. For added convenience, the sagittal saw head can be rotated in 45° increments to provide additional flexibility during cases. The MicroFree® System also features the modular Mini-Driver™ for pin/wire-driving and drilling as well as the system's latest battery-powered, pencil-grip wire driver handpiece.

SYSTEM INCLUDES:

- 6 HANDPIECES
- MICROFREE® LITHIUM POWER UNIT
- MINI-DRIVER
- SMALL LITHIUM BATTERY
- CUSTOMIZABLE STERILIZATION TRAYS


HALL®
MICROFREE®
SMALL BONE
POWER SYSTEM



Hall® Mini-Driver



Small Lithium Battery

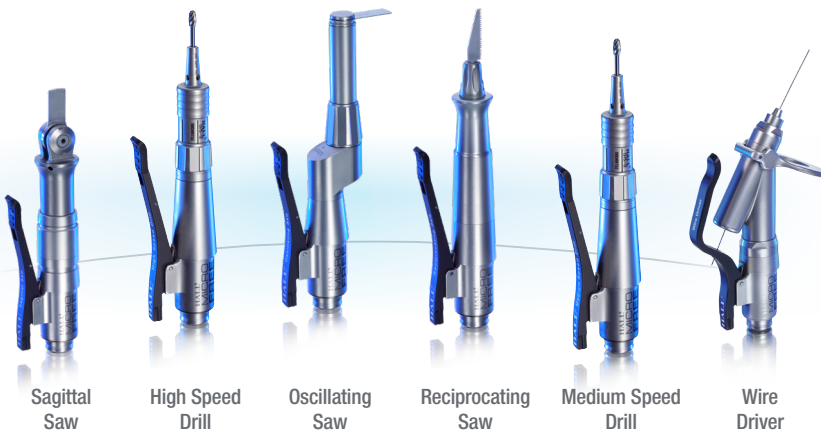

Sagittal Saw

High Speed Drill

Oscillating Saw

Reciprocating Saw

Medium Speed Drill

Wire Driver

Powerful High-Performance Lithium-Ion Technology

A cordless system is only an asset if it is powerful enough to get the job done. Our high-performance lithium-ion technology provides the torque, power, and performance needed.




Lithium Power Unit

HALL®
MICROFREE®
SMALL BONE
POWER SYSTEM

Facility Benefits: One System That is Easy to Maintain and Upgrade

Autoclavable for Efficient and Cost-Effective Sterile Processing

Both the UL-approved batteries and handpieces of the MicroFree® System are autoclavable, eliminating the need for more expensive STERRAD battery sterilization. This efficient system allows the batteries, handpieces, and attachments to be sterilized together in the same tray at the same time – simplifying the process and saving OR staff time and energy.

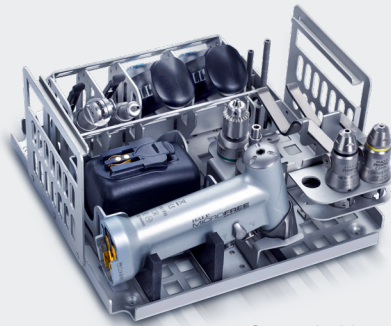

Customizable
Sterilization Trays



MicroFree
Handpiece

IPX
6

*Protected against high pressure
water stream from any angle.*

MicroFree
Handpiece

IPX
8

IPX
9

*Protected against continuous submersion
in water over one meter deep.*

IPX 6, IPX 8, and IPX 9 - Protection from Moisture

Just like our large bone Hall® Titan® handpiece, the MicroFree handpiece is rated IPX 6, IPX 8, and IPX 9 to help protect the system against water damage. The UL-approved hermetically sealed power units are fully autoclavable, allowing the whole system to be sterilized together.

Safety – Designed to Reduce the Risk of Infection

By eliminating unsecured power cords from compromising the sterile field, the MicroFree® System may reduce those potential risks that can contribute to surgical site infections (SSIs).



Combined with the Hall® Titan® System for Large Bone, the MicroFree® System can outfit your facility with a power solution that is completely autoclavable.

Backward and Forward Compatibility for Easy Upgrades

Often when other companies roll out new systems, they are intentionally incompatible with their old technology – forcing expensive upgrades when your existing equipment is still working. Like most Hall® Systems, MicroFree® is designed to be both backward and forward compatible – delivering greater long-term value for your investment.

High Quality, Flexible Service Options

As a FDA-regulated and ISO 13485 certified facility, CONMED service provides fast, cost-effective, and high-quality service for every component of the MicroFree® System. We also offer complete flexibility in our service solutions and multiple program levels that can be tailored to meet the specific needs of your organization's surgical schedule, inventory, and budgetary requirements.

Experience the Freedom of MicroFree® – Schedule an Evaluation Today

To take the next step toward discovering the comfort and convenience of cordless small bone power, schedule an evaluation by calling **1-800-237-0169** to locate your local CONMED Sales Representative.

For more information, visit www.CONMED.com/MicroFree.

HALL®
MICROFREE®
SMALL BONE
POWER SYSTEM

Ordering Information

Description	Catalog Number
Hall® MicroFree® Handpieces	
Hall MicroFree High Speed Drill	PR08000SB
Hall MicroFree Medium Speed Drill	PR08100SB
Hall MicroFree Sagittal Saw	PR08200SB
Hall MicroFree Oscillating Saw	PR08300SB
Hall MicroFree Reciprocating Saw	PR08400SB
Hall MicroFree Pencil-Grip Wire Driver	PR08600SB
Hall® MicroFree® Lithium Power Unit	
Hall MicroFree Lithium Power Unit	L3500SB
Hall® MicroFree® Mini-Driver	
Hall MicroFree Mini-Driver	PR08500SB
Hall® Small Lithium Battery	
Hall Small Lithium Battery	L3000SM
Hall® MicroFree® Adapter and Charger	
Hall MicroFree Power Unit Adapter	L3500
Hall Lithium Battery Charger	L3000
Hall® MicroFree® Sterilization Containers	
1/2 Size Small Bone Insert	PR08003T
1/2 Size Procedure Specific Insert	PR08001T
1/2 Size Rigid Container	TR12R
3/4 Size Rigid Container	TR34R
Accessories: Bur Guards	
Medium Bur Guard	1375-012
Long Bur Guard	1375-011
Extra-Long Bur Guard	1375-023